

BD™ Blunt Fill and BD™ Blunt Filter Needles


Quick reference guide



1. Attachment

Attach the BD™ Blunt Fill or BD™ Blunt Filter needle to any standard luer lock or luer slip hypodermic syringe.

Twist until firmly seated to connect to a luer lock syringe. Pull shield straight off the needle to avoid damaging the needle point.

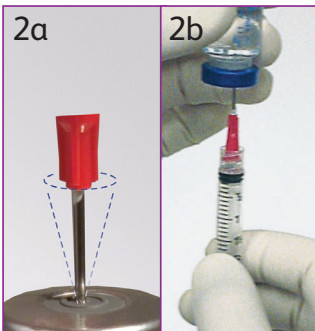


2. Aspiration

2a. Puncture stopper with the BD™ Blunt Fill needle.

Note: For glass ampules, use the BD™ Blunt Filter needle.

2b. Draw up medication in accordance with your institution's protocol.

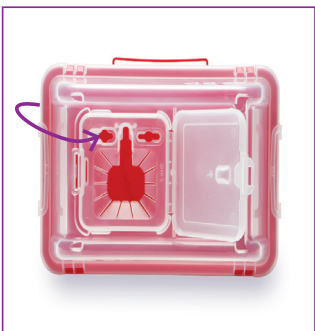

3. Disposal

Grip the needle hub using the v notch available on approved sharps container and immediately discard.

Remember

- The BD™ Blunt Fill and BD™ Blunt Filter needles are not to be used for skin injections.
- Follow your institution's protocol for aseptic technique.
- Follow your institution's policy for safe disposal of all medical waste.
- For glass ampules, use the BD™ Blunt Filter needle.
- The BD™ Blunt Filter needle is a one-way filter and should be removed once fluid is filtered into the syringe.
- To connect the BD™ Blunt Fill or BD™ Blunt Filter needle to a luer slip syringe, push firmly.

BD™ Blunt Fill Needle



BD™ Blunt Filter Needle



Note hub color difference between blunt filter vs. blunt fill needle